

Manufactured by:

ResMed Ltd | Elizabeth MacArthur Drive Bella Vista NSW 2153 Australia
Tel: +61 (2) 8884 1000 or 1 800 658 189 (toll free) Fax: +61 (2) 8883 3114 Email: reception@resmed.com.au

ResMed Offices:

ResMed Corp 14040 Danielson Street Poway CA 92064-6857 USA
Tel: +1 (858) 746-2400 or 1-800-424-0737 (toll free) Fax: +1 (858) 746-2900 Email: reception@resmed.com

ResMed (UK) Ltd 65 Milton Park Abingdon Oxfordshire OX14 4RX UK
Tel: +44 (1235) 862 997 Fax: +44 (1235) 831 336 Email: reception@resmed.co.uk

ResMed Asia Pacific Ltd 97 Waterloo Road North Ryde NSW 2113 Australia
Tel: +61 (2) 9886 5000 or 1 800 991 900 (toll free) Fax: +61 (2) 9889 1471 Email: reception@resmed.com.au

ResMed Finland Niittykatu 6 FIN 02200 Espoo Suomi
Tel: +358 9 8676820 Fax: +358 9 86768222 Email: reception@resmed.fi

ResMed GmbH & Co. KG Rudolfstraße 10 D-41068 Mönchengladbach Deutschland
Tel: +49 (0) 2161-3521-0 (Reception), +49 (0) 180 22 22 668 (Service-Telefon; 0.06 s/Anruf),
+49 (0) 180 22 66 888 (Wartungstelefon; 0.06 s/Anruf) Fax: +49 (0) 2161-3521-1499 Email: reception@resmed.de

ResMed Hong Kong Ltd Room 1714 Miramar Tower 132-134 Nathan Road Tsim Sha Tsui Hong Kong
Tel: +852 2366 0707 Fax: +852 2366 4546 Email: reception@resmed.com.hk

ResMed Japan Nihonbashi Hisamatsu Bldg. 4F 2-28-1 Nihonbashi-Hamacho Chuo-Ku, Tokyo 103-0007 Japan
Tel: +81 (3) 3662 5056 Fax: +81 (3) 3662 5040

ResMed Malaysia Sdn Bhd Suite E-10-20, Plaza Mon't Kiara No. 2, Jalan 1/70C Mon't Kiara 50480 Kuala Lumpur
Malaysia Tel: +60 3 6201 7177 Fax: +60 3 6201 2177 Email: reception@resmed.com.my

ResMed NZ Ltd PO Box 51-048 Pakuranga Auckland New Zealand
Tel: +0800 737 633 (NZ toll free) or +64 274 737 633 Fax: +0800 737 634 (NZ toll free) or +64 9 239 0193
Email: reception@resmed.co.nz

ResMed SA Parc de la Bandonnière 2 rue Maurice Audibert 69800 Saint-Priest France
Tel: +33 (0) 4 37 251 251 Fax: +33 (0) 4 37 251 260 Email: reception@resmed.fr

ResMed Singapore Pte Ltd 238A Thomson Road # 12-03/04 Novena Square Tower A Singapore 307684
Tel: +65 6284 7177 Fax: +65 6284 7787 Email: reception@resmed.com.sg

ResMed Spain SL C/ Arturo Soria, 245 28033 Madrid España
Tel: +34 (93) 590 8154 Fax: +34 (93) 590 8153 Email: angelo@resmed.es

ResMed Sweden AB Industrigatan 2 S-461 37 Trollhättan Sverige
Tel: +46 520 420 110 Fax: +46 520 397 15 Email: reception@resmed.se

Labhardt AG Thannerstrasse 57 CH-4054 Basel Schweiz
Tel: +41 (061) 307 9711 Fax: +41 (061) 307 9722 Email: info@labhardt.ch

US DESIGNATED AGENT: ResMed Corp
EU AUTHORISED REPRESENTATIVE: ResMed (UK) Ltd

Mirage Kidsta Nasal Mask
Protected by patents: AU 710733, AU 712236, AU 724360, DE 29724224, US 6112746, US 6357441, US 6439230, US 6561190, US 6561191, US 6581602, US 6634358. Other patents pending.
Protected by design registrations: AU 155372, CN 3302716.1, DE 40300897, FR 03 0845, GB 3010608, GB 3014182, GB 3014183, GB 3014184, GB 3014185, GB 3014186, GB 3014187, GB 3014188, GB 3014189, JP 1200870, JP 1201246, NZ 403295, SE 77095, US D485905. Others pending.

Mirage, Kidsta and Mirage Kidsta are trademarks of ResMed Ltd.
Dove is a registered trademark of Unilever.

© 2004 ResMed Ltd.

RESMED

Mirage Kidsta™ Nasal Mask

Reorder No. 61809/2 04 08



USER'S GUIDE
English



www.resmed.com

Waking people up to sleep


© 2004 ResMed Ltd.

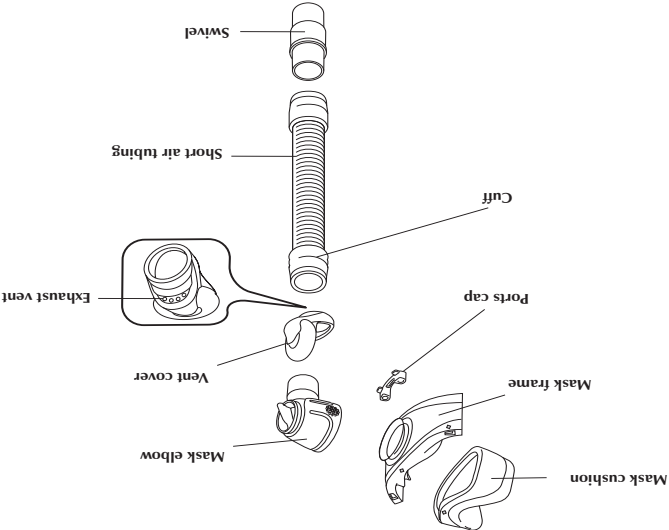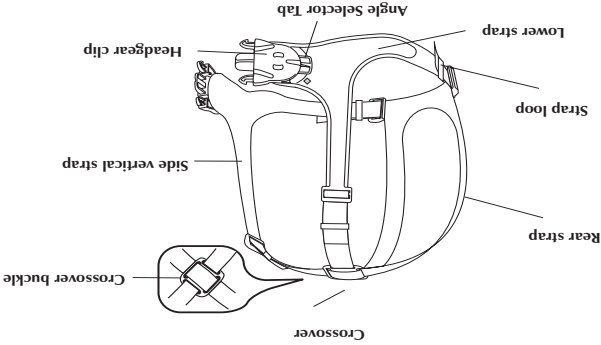

Mask Parts

B

© 2004 ResMed Ltd.






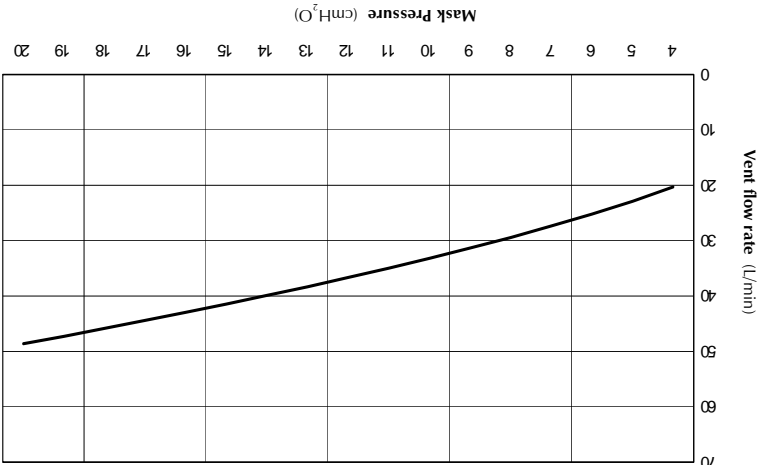


Pressure-Flow Curve

A


Note: For complete instructions, read the text section in conjunction with the illustrations on this sheet.

C Fitting

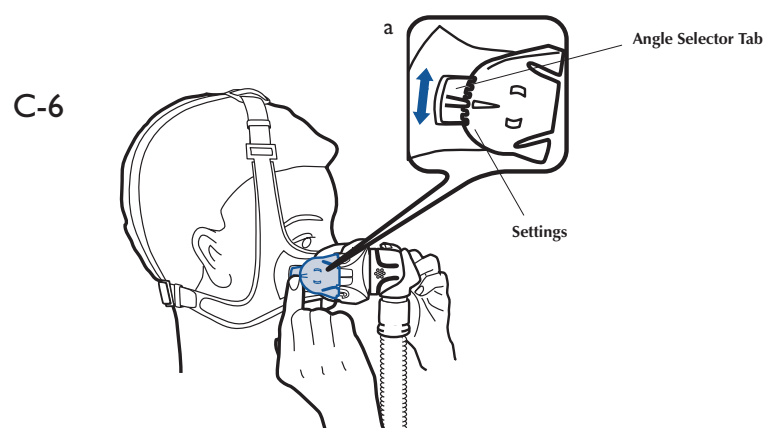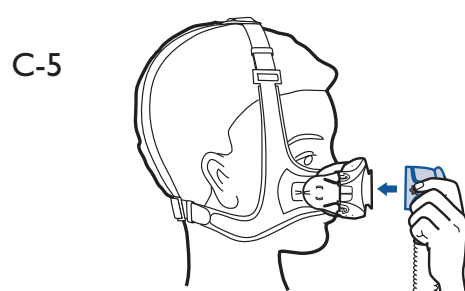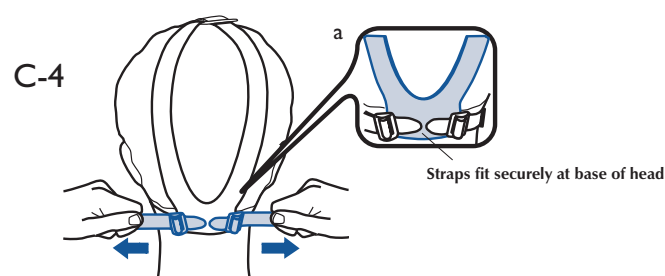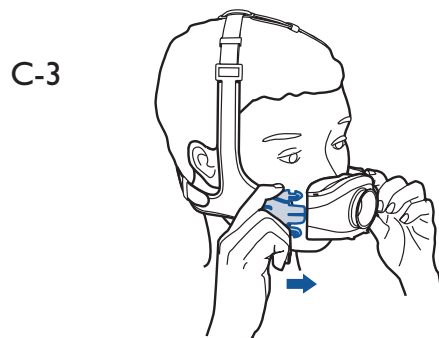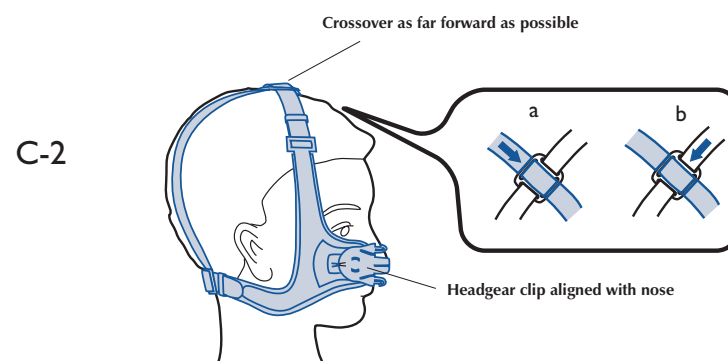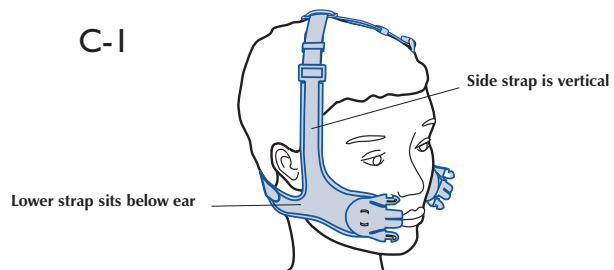


© 2004 ResMed Ltd.

D Disassembling

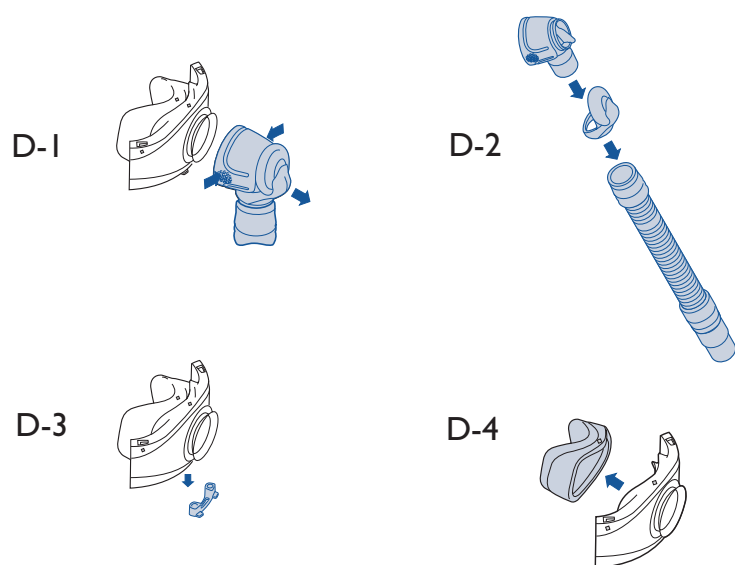


© 2004 ResMed Ltd.

E Reassembling

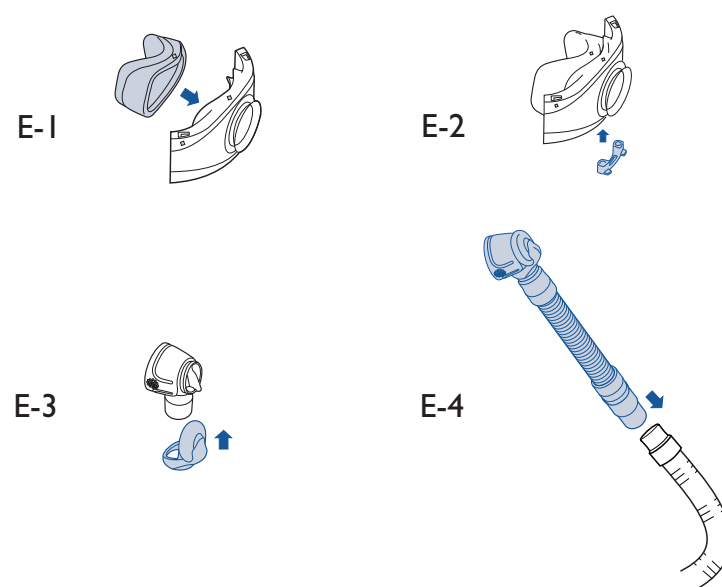

© 2004 ResMed Ltd.

Note: The headgear does not need to be disassembled for cleaning. Headgear reassembly instructions are provided here as additional information.

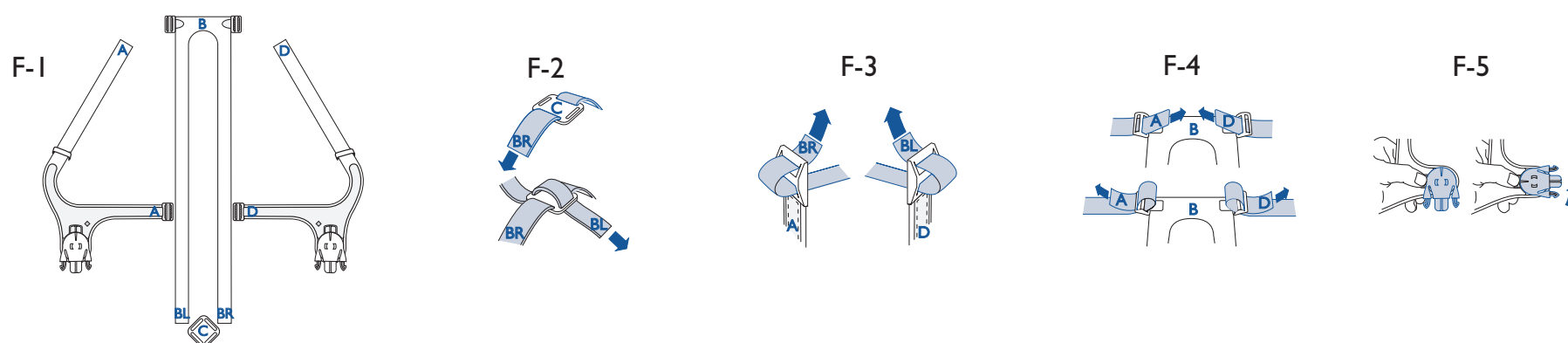

© 2004 ResMed Ltd.

Mirage Kidsta™ Nasal Mask

The **MIRAGE KIDSTA™ NASAL MASK** is intended for multipatient use for children aged seven and over, prescribed continuous positive airway pressure (CPAP) or bilevel therapy in hospital, clinic and/or home environments.

MEDICAL INFORMATION

Note: *The mask system does not contain latex materials. If you have ANY reaction to any part of the mask system, discontinue use and consult your sleep therapist.*



WARNINGS AND CAUTIONS

- Your physician should determine the degree of arousability before commencing CPAP treatment. Use of secondary sensing devices such as apnoea monitors and transcutaneous carbon dioxide monitors is advisable. These will help determine whether nasal CPAP treatment should be applied.
- Ensure that the nasal cushion does not block the nostrils at any time during fitting or use.
- The **MIRAGE KIDSTA** is designed to provide a stable, effective air seal with minimum headgear tension. High headgear tension may lead to facial deformities with long-term use. Regular clinician review of growth patterns and

the use of different styles of mask is recommended.

- The mask contains small parts which could result in a choking risk. Your clinician should assess the suitability of the selected mask.
- This mask should be used only with flow generators recommended by your physician or respiratory therapist. A mask should not be used unless the flow generator is turned on and operating properly. The exhaust vent should never be blocked.

Explanation of the Warning: Flow generators are intended to be used with special masks with connectors which may have vent holes or separate exhaust vents to allow continuous flow of air out of the mask. When the flow generator is turned on and functioning properly, new air from the flow generator flushes the exhaled air out through the attached mask exhaust vent. However, when the flow generator is not operating, enough fresh air will not be provided through the mask, and exhaled air may be rebreathed. Rebreathing of exhaled air for longer than several minutes can in some circumstances lead to suffocation. This warning applies to most models of flow generators.

- At low CPAP pressures, the flow through the exhaust vent may be inadequate to clear all exhaled gas from the tubing. Some rebreathing may occur.

Explanation of the Warning: The exhaust vent flow is lower at lower mask pressures (see Pressure/Flow curve—Figure A). To avoid rebreathing it is recommended that a mask pressure of 4 cmH₂O or more be used.

- The **MIRAGE KIDSTA** features an elbow that can be easily removed while the mask is fitted. The mask should not be used with the elbow attached unless the flow generator is turned on and operating properly.
- If oxygen is used with the flow generator, the oxygen flow must be turned off when the flow generator is not operating.

Explanation of the Warning: When the flow generator is not in operation, and the oxygen flow is left on, oxygen delivered into the ventilator tubing may accumulate within the flow generator enclosure. Oxygen accumulated in the flow generator will create a risk of fire. This warning applies to most types of flow generators.

- Please refrain from smoking while oxygen is in use.

Note: *At a fixed rate of supplemental oxygen flow, the inhaled oxygen concentration will vary, depending on the pressure settings, patient breathing pattern, mask selection and leak rate.*

PARTS OF THE MASK

See **Figure B**.

FITTING THE MASK

See **Figure C**.

Note: *Before fitting the mask, refer to the section “Cleaning the Mask in the Home”.*

- 1 Place the headgear on the head. Adjust the side vertical straps so that the lower strap sits just below the ear (not touching it) and at the base of the head (**Fig. C-1**).
- 2 The headgear clip should be parallel with the nose and not in the line of vision. Check that the crossover is as far forward as possible (**Fig. C-2**). To adjust the crossover to a forward position, slide the crossover buckle along one strap (**Fig. C-2a**) and then along the other strap (**Fig. C-2b**). Alternate the movement on each side until the crossover is as far forward as possible.
- 3 Ensure that the crossover and the rear straps lie flat against the head.
- 4 Place the mask on the nose (**Fig. C-3**). Clip the mask to the headgear.
- 5 Adjust the lower straps by pulling firmly (**Fig. C-4**) so that the rear straps fit securely around the base of the head (**Fig. C-4a**).

Note: *Ensure that the nostrils are not blocked by the cushion.*

- 6 Connect one end of the main air tubing to your flow generator and the other end to the swivel.

Note: *The swivel, short air tubing and mask elbow are supplied assembled.*

- 7 Turn your flow generator on.
- 8 Connect the mask elbow to the mask frame. To do this, hold the sides of the mask elbow (**Fig. C-5**) and firmly push the elbow onto the mask frame until it clicks into place.

Check that you have connected the mask elbow correctly by rotating it. Be careful not to squeeze the sides of the elbow. If the mask elbow is fitted correctly, it should rotate freely and not separate from the mask frame.

- 9 Adjust the angle of the mask by selecting the setting that provides a comfortable fit with minimum or no leak. To test each setting, hold down the Angle Selector Tab (**Fig. C-6**) and move the headgear clip up or down (**Fig. C-6a**) as required. You may need to raise the headgear clip when moving it from one setting to the next. The mask can be tilted at five different positions. If you feel excessive pressure on the middle of the nose, tilt the mask down. Alternatively, if you feel excessive pressure on your upper lip, tilt the mask up.
- 10 If you feel air leaks around the cushion, reposition the mask on your face to achieve a better seal. Take time to adjust all straps and the position of the mask so that you can achieve a secure and comfortable fit.

Notes:

- If you are unable to resolve any leaks you may need to try another cushion size.
- Excess headgear length may be trimmed to a suitable length if required. Leave at least 30 mm free length on the straps.
- Select "Standard" as the mask option when using the MIRAGE KIDSTA with ResMed flow generators that have mask setting options.

REMOVING THE MASK

To remove the MIRAGE KIDSTA, squeeze one of the headgear clips. This releases one side of the headgear from the mask. Pull the mask and headgear away from your head.

Notes:

- Do not pull the headgear off your head without releasing the headgear clips first.
- You can detach the elbow from the mask frame while wearing the mask. This is useful when therapy needs to be stopped for a short time; for example if you need to go to the bathroom.
- Check that you can remove the headgear clip, the quick release swivel or the headgear unaided.

DISASSEMBLING FOR CLEANING

See **Figure D**.

- 1 Detach the mask elbow. To do this, squeeze the sides of the mask elbow at the finger grips while pulling it away from the mask frame (**Fig. D-1**).
- 2 Disconnect and separate the main air tubing, short air tubing, mask elbow and vent cover (**Fig. D-2**).

Notes:

- To prevent damage to the short air tubing, disconnect it by holding it at the cuff.
 - The swivel does not need to be disconnected from the short air tubing.
- 3 Detach the headgear from the mask frame by squeezing and releasing both headgear clips.
 - 4 Remove the ports cap from the mask frame (**Fig. D-3**).
 - 5 Firmly pull the mask cushion away from the mask frame (**Fig. D-4**).

CLEANING THE MASK IN THE HOME

Notes:

- Wash the headgear before first time use, as the dye may run. It is important to always handwash the headgear.
- The headgear may be washed without being disassembled.

DAILY/AFTER EACH USE

Handwash the mask components in lukewarm, soapy water. Use pure soap (an example of pure soap is Dove® soap, unscented; made by Unilever). Rinse the components well and allow them to air dry out of direct sunlight.

WEEKLY

Handwash the headgear in warm, soapy water. Use pure soap. Rinse well and allow it to air dry out of direct sunlight.



CAUTIONS

- Do not use solutions containing bleach, chlorine, alcohol, aromatics, moisturisers, antibacterial agents, or scented oils to clean any part of the mask system or air tubing. These solutions may cause damage and reduce the life of the product.
- Do not expose any part of the mask system or tubing to direct sunlight as it may deteriorate.
- If any visible deterioration of a mask component is apparent (cracking, crazing, tears etc), the mask component should be discarded and replaced.

CLEANING THE MASK BETWEEN PATIENTS

The mask system (excluding the swivel and short tube) is validated for multipatient use. Clinicians should refer to the *Masks Disinfection/Sterilisation Guide* for details.

To clean the headgear between patients, use the method described in the previous section “Cleaning the Mask in the Home”.

REASSEMBLING THE MASK

See **Figure E**.

- 1 Insert the mask cushion into the groove of the mask frame (**Fig. E-1**). Fit the top section first and then work your way around, pushing the mask cushion into the mask frame.

Note: The diamond markers on the mask frame and on the mask cushion (one on top, two at the bottom) must line up.

- 2 Push the ports cap onto the mask ports (**Fig. E-2**).
- 3 Attach the vent cover to the mask elbow. To do this, insert the elbow into the ring of the vent cover (**Fig. E-3**). Push the ring to the top so that the vent cover is over the back of the elbow. Then, press the lip of the vent cover to ensure that it is securely seated in the groove of the elbow.
- 4 Connect the short air tubing to the mask elbow.
- 5 Connect the swivel to the main air tubing, ready for next use (**Fig. E-4**).

ADDITIONAL INFORMATION

REASSEMBLING THE HEADGEAR

Note: *You do not need to disassemble the headgear for cleaning. However, if you choose to do so, the reassembly procedure is explained in this section.*

When disassembled, the headgear has four parts:

A: Left strap, B: Middle strap,
C: Crossover buckle, D: Right strap

- 1 Lay out the four parts of the headgear as shown in **Figure F-1**.
- 2 As shown in **Figure F-2**, thread the straps BR and BL through the crossover buckle (Part C).
- 3 Thread strap BR through the buckle on Part A. Thread strap BL through the buckle on Part D (**Fig. F-3**).
- 4 Thread the straps of Part A and Part D through the buckles on Part B (**Fig. F-4**).
- 5 With the headgear clip facing down, fit the headgear clip onto the hole on the headgear. While pressing the headgear tab down, rotate the headgear clip away from the tab to lock it into place (**Fig. F-5**). Repeat this with the other headgear clip.

TECHNICAL SPECIFICATIONS

PRESSURE–FLOW CURVE

See **Figure A**.

DEAD SPACE INFORMATION

Dead space is the empty volume of the mask up to the swivel.

The dead space of the mask for all cushion sizes is 81 mL.

REPLACEMENT PARTS

Part No.	Item
61020	Headgear Assembly, Small
60013	Cushion, Small
61025	Headgear clip - blue (2 pk.)
61026	Headgear clip - blue (10 pk.)
60023	Crossover buckle - blue
60024	Crossover buckle - blue (10 pk.)
60924	Mask Frame
60925	Elbow Assembly
60014	Swivel Assembly (10 pk.)
60926	Inlet Tube and Swivel Assembly
60927	Inlet Tube and Swivel Assembly (10 pk.)
16570	Ports cap (2 pk.)
16571	Ports cap (10 pk.)
61011	Mask System, Small

LIMITED WARRANTY

ResMed warrants that your ResMed mask system, including its components: mask frame, cushion, headgear, tubing and other accessories ("Product") shall be free from defects in material and workmanship for a period of ninety (90) days from the date of purchase by the initial consumer.

For further details refer to the Warranty.

Note: *Translations are based on the English version in this guide. However, there may be minor variations based on the disinfection/sterilisation procedures applicable in some countries. ResMed has reviewed these changes and all procedures included here have been validated.*



Component of: 61809/2 04 08